

**Never remove  
an inflated  
balloon again**



NEW  
INNOVATIVE  
DESIGN

Try the next generation of Foley catheters

ORDER NOW



**One small step  
for vets. One giant  
leap for animals**

NEW  
INNOVATIVE  
DESIGN



With this innovative, easy-to-use solution, you'll never have to say "oh no!" again. Let your patients wake up without the agony of a catheter removal gone wrong.



Visit [kruise.com/safvetfoley](https://kruise.com/safvetfoley)  
and learn more about the new  
KRUUSE SafeVet® Foley Catheter

kruise.com



Gentle Treatment

KRUUSE 08/2024 - 21257/INT

**THE  
WORST  
MISTAKES**



**...are the ones that  
could have been  
easily avoided**



Gentle Treatment



# "Did I empty it?"

You're not alone: Many vets have tried to remove a Foley catheter with the balloon still inflated. We're all busy. Everyone can forget. But now, no one has to.

The new KRUUSE SafeVet® Foley Catheter eliminates all doubt, so you can always remove it safely.



See how it works

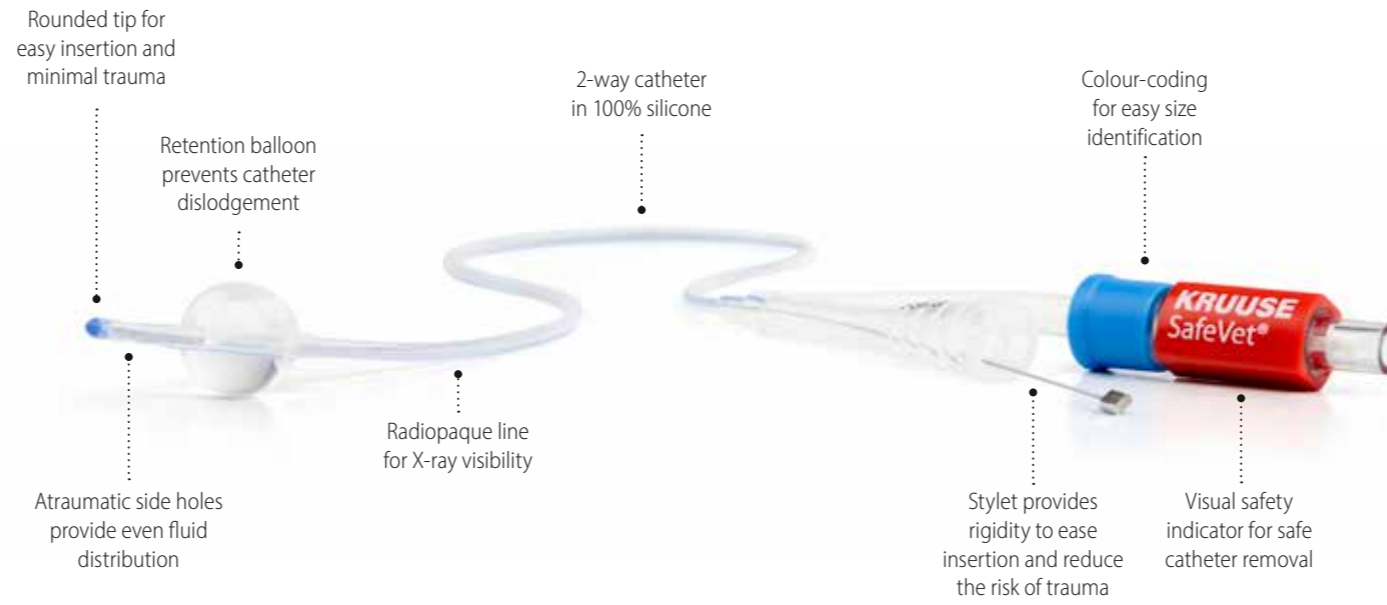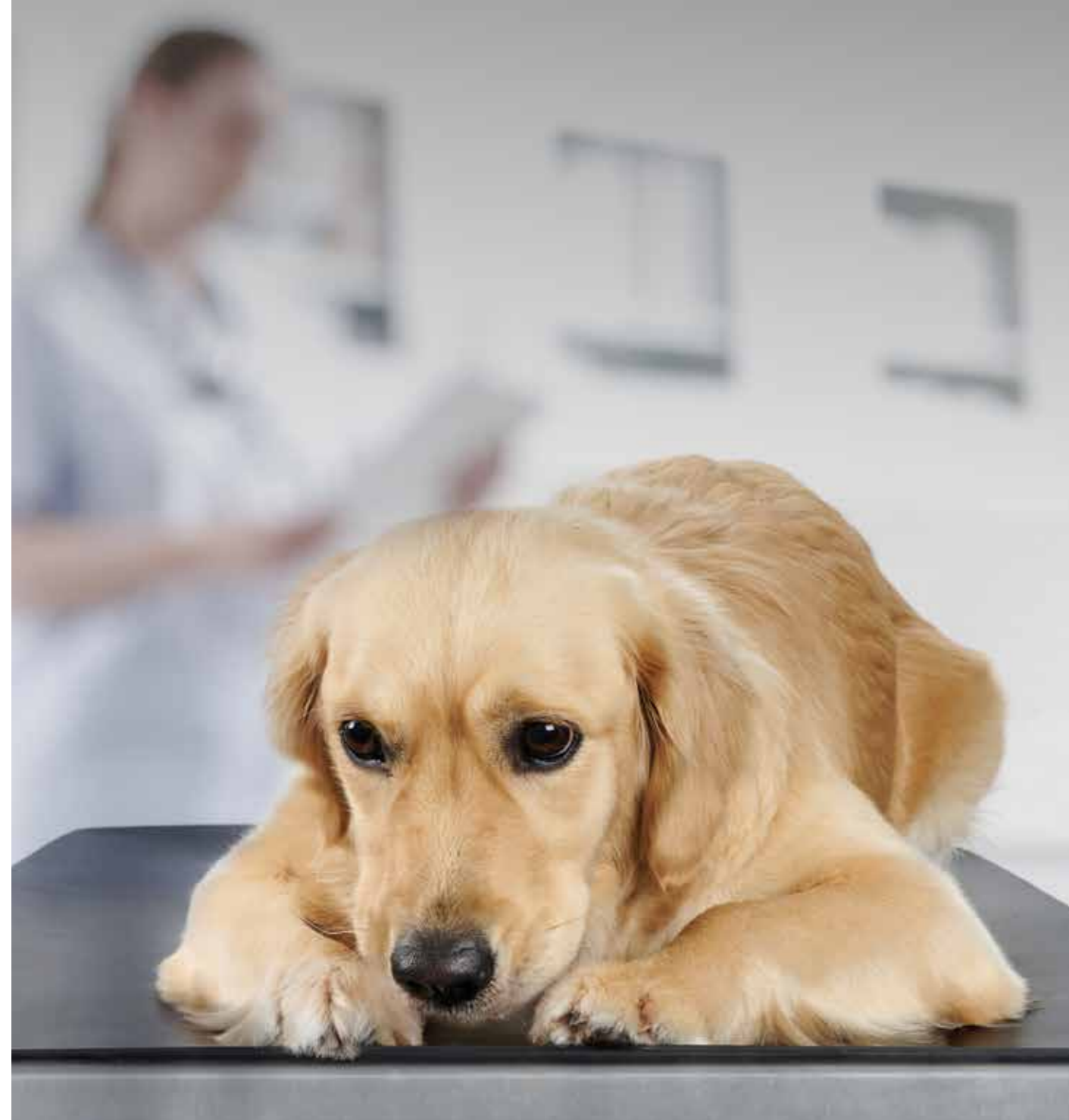

## Eliminate doubt. Eliminate mistakes

### A choice for compassionate care

The KRUUSE SafeVet® valve serves as a safety mechanism for removal of the KRUUSE SafeVet® Foley Catheter. To initiate catheter removal, simply unscrew the KRUUSE SafeVet® valve. Wait until the catheter is fully emptied for safe removal.

Choosing Foley catheters for long-term catheterisation demonstrates a veterinary practice's commitment to providing compassionate, state-of-the-art care.

KRUUSE SafeVet® Foley Catheters are made from highly biocompatible silicone, ensuring maximum patient comfort.

### Designed for enhanced safety and comfort

- Reduces unnecessary pain
- Eliminates doubt
- Prevents mistakes



### KRUUSE SafeVet® Foley Catheter

Cat. No	Type	Sterile	Colour	Amount cc (ml) in retention balloon	Size	Package
273832	Side holes, with wire	Yes	● Green	3 cc (ml)	2.0 mm x 30 cm / 6 Fr x 12"	5/pk
273833	Side holes, with wire	Yes	● Green	3 cc (ml)	2.0 mm x 55 cm / 6 Fr x 22"	5/pk
273834	Side holes, with wire	Yes	● Blue	3-5 cc (ml)	2.7 mm x 30 cm / 8 Fr x 12"	5/pk
273835	Side holes, with wire	Yes	● Blue	3-5 cc (ml)	2.7 mm x 55 cm / 8 Fr x 22"	5/pk
273830	Side holes, with wire	Yes	● Blue	3-5 cc (ml)	2.7 mm x 70 cm / 8 Fr x 28"	5/pk
273836	Side holes, with wire	Yes	● Black	3-5 cc (ml)	3.3 mm x 30 cm / 10 Fr x 12"	5/pk
273837	Side holes, with wire	Yes	● Black	3-5 cc (ml)	3.3 mm x 55 cm / 10 Fr x 22"	5/pk
273831	Side holes, with wire	Yes	● Black	3-5 cc (ml)	3.3 mm x 70 cm / 10 Fr x 28"	5/pk

## Just like you're used to. Without the fear of inflicting pain

You don't have to change the way you do things. All you have to do is look at the valve.



### Try the new multi-purpose adaptor

KRUUSE Foley Adaptor with Closing Cap can be used to seal the catheter drainage funnel to prevent leakage or as a connection to a urine bag or syringe to allow flow.



Cat. No	Type	Sterile	Colour	Package
273838	KRUUSE Foley Adaptor with Closing Cap	Yes	● Light blue	5/pk